



Alere™ Toxicology Services

OrAlert™ Oral Fluid Drug Test Device


In addition to offering Point of Collection Testing (POCT) devices that detect drugs in urine, Alere can supply its clients the means to detect drugs in oral fluid.

The OrAlert™ Oral Fluid Drug Test Device is a method for the noninvasive, instant detection of multiple drugs in human oral fluid.

Using a multi-drug, multi-line lateral flow immunoassay, the OrAlert™ enables the qualitative detection of drugs and drug metabolites. The device's use of oral fluid eliminates the need for special equipment and same-gender collections. Additionally, because oral fluid is difficult to alter or adulterate, the specimen's integrity can be assured.

Configuration for Product Number DSF-765-011:

- THC – Marijuana
- COC – Cocaine
- OPI – Opiates, including morphine and codeine
- AMP/MAMP – Amphetamines/Methamphetamines
- PCP – Phencyclidine

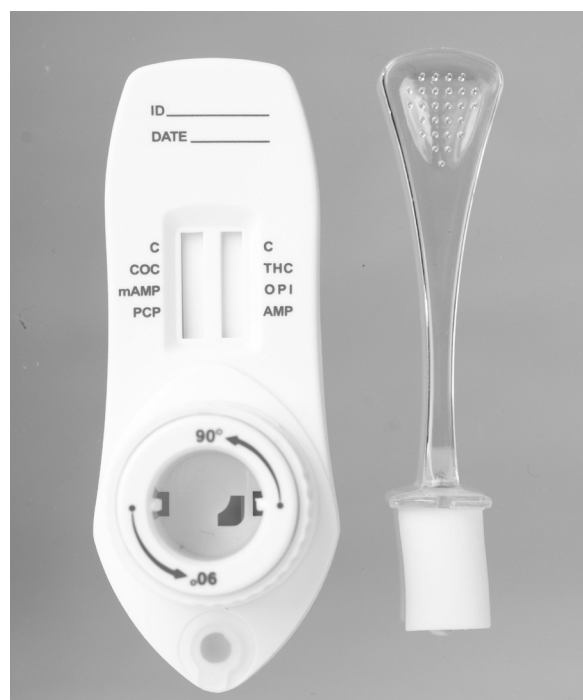

Features:

- No additional special equipment or environment is required for use of the device.
- Donors are not required to provide specimens in privacy, so collections may be performed by the opposing gender.
- Test devices are self-contained and portable and can be sent directly to our laboratory for confirmation.
- Test controls are built into the device.
- Results are available in 5 minutes.
- Results appear clear and are easy-to-read.
- The devices may be stored at room temperature for 24 months.
- Devices must be ordered in increments of 25.
- Devices are for forensic use only.